LIDÉN BASIC

LEP-185

Version 1.0 2026-06-04

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### Introduction

---

**For the best performance and weighing results, take care in the daily use and maintenance of your scale.**

**For long durability and best results ensure that:**

- **DO NOT** overload the scale. This will damage the load cell and void the warranty.
- **DO NOT** leave the scale with a load on it for an extended period. This reduces the accuracy of the scale and may shorten the life of the load cell.
- **DO NOT** try to repair the scale yourself, in case of problems, contact your local representative.
- Avoid exposing the scale to temperature differences. Do not place the scale in direct sunlight or near cold areas.
- Make sure that the scale is on a stable surface and that it is not subject to vibrations. If your scale is equipped with adjustable feet and spirit levels, adjust if necessary.
- Do not place the scale near large electrical appliances, such as welders or large motors. Do not leave batteries in the scale if it is not in use for a long time.
- Avoid placing the scale in a humid environment where condensation can occur. Avoid exposing the scale to direct contact with water.
- **DO NOT** place the scale near windows, doors, or fans where drafts can affect the weighing result.
- Keep the scale clean. Wipe, if necessary, with a damp cloth.
- Be sure to use the original adapter/charger and battery if the scale is equipped with it.

### **WARNING – Safety Information**

---

**Read this information carefully before use. Incorrect use may cause personal injury, inaccurate measurement results, or damage to the product.**


- Do not use the body analysis function if you have a pacemaker or other implanted medical device.
- During pregnancy, use the scale for weight measurement only.
- Do not use the scale with wet feet or on a wet weighing platform.
- Never load the scale with more than 180 kg.
- Never jump on the scale or use it as a step stool or support.
- Do not allow children to play with the scale. Keep packaging materials out of reach of children.
- Stop using the scale immediately if it emits smoke, unusual odors, or excessive heat.
- Avoid contact with leaking battery fluid. In case of contact with skin or eyes, rinse immediately with plenty of water.



### Use

---

#### Before use

Remove the plastic strip from the battery. Battery compartment (7) under the scale. The scale is now ready to use.

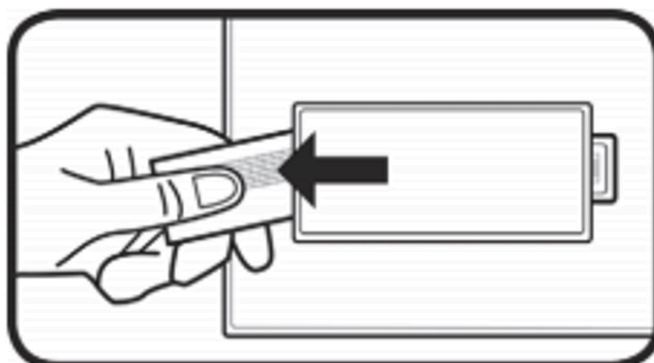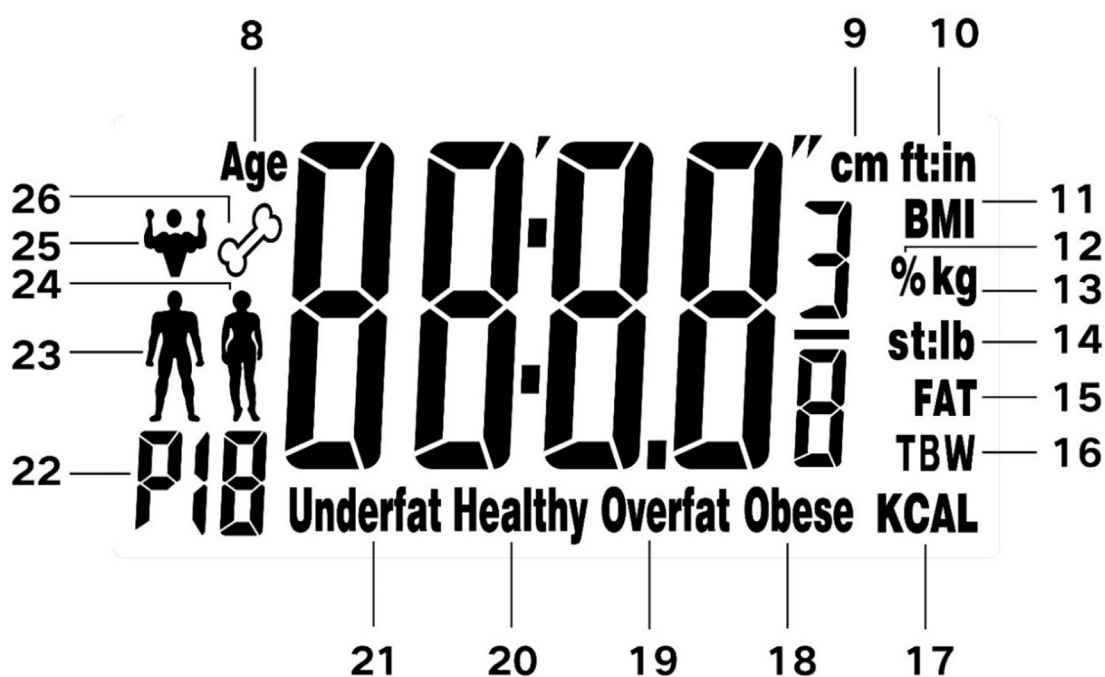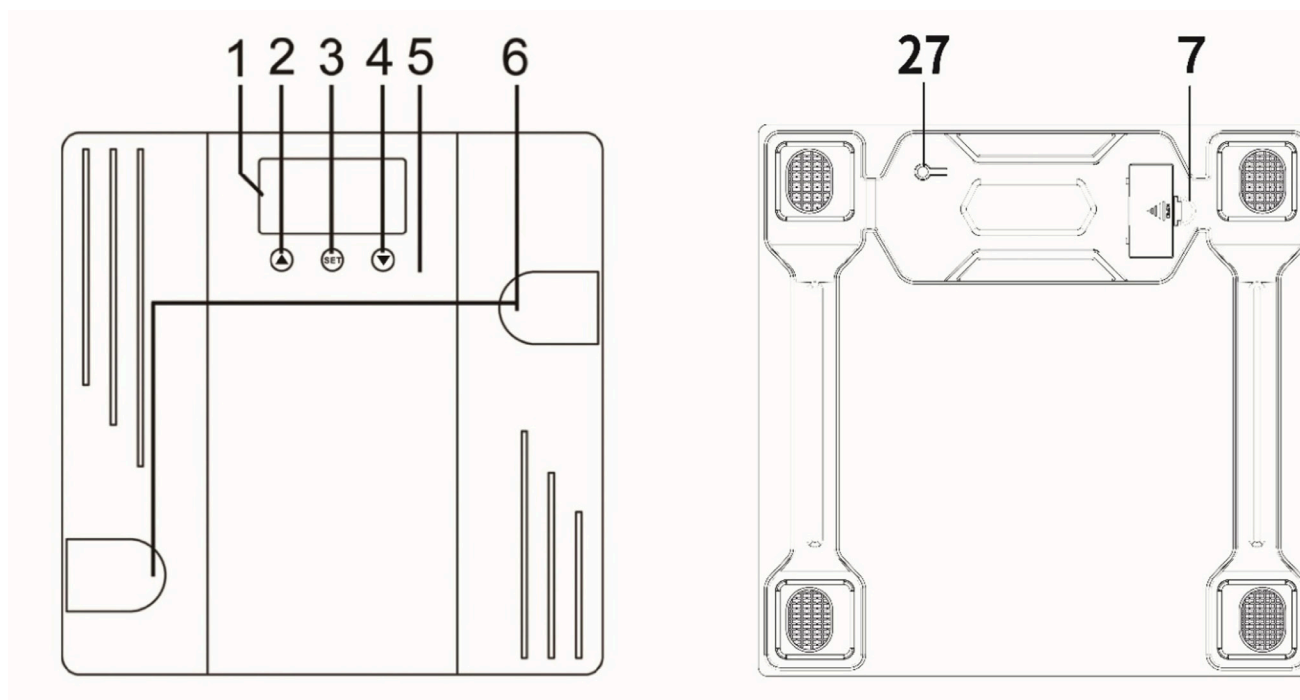- |                                |   |                                 |
|--------------------------------|---|---------------------------------|
| 1. LCD Display                 | 11. BMI                                       | 20. Normal weight               |
| 2. Button ( <b>UP</b> )        | 12. Percent value (of bodyfat, water, muscle) | 21. Underweight                 |
| 3. <b>SET</b> button           | 13. Weighing unit kg                          | 22. Memory for 12 user profiles |
| 4. Button ( <b>DOWN</b> )      | 14. St:lb unit = stones:pounds.               | 23. Man                         |
| 5. Scale's platform            | 15. Bodyfat in percent                        | 24. Woman                       |
| 6. Elektrodes                  | 16. Water amount in percent                   | 25. Muscle mass in percent      |
| 7. Battery compartment         | 17. Recomend calorie intake / day             | 26. Bone mass                   |
| 8. Age indication              | 18. Obesity                                   | 27. Unit                        |
| 9. Body length in cm           | 19. Overweight                                |                                 |
| 10. Body length in feet/inches |   |                                 |

## Start

Due to the 'Step on power' function, the digital scale turns on automatically when you step on the scale plate (5) and it starts measuring your weight immediately.

If **0.0** is not displayed at startup, wait about three seconds for reset.

## Change unit

In its original state, the scale is set to kg (13) as the unit of weight. If you want to change the weight unit, press button (27) on the back. The units change in turn after each press.

## Weighing

1. Stand on the scale (5) so that your weight is evenly distributed. Stay as still as possible during the weighing process, otherwise the measurement result may be misleading.
2. The measurement result is shown on the LCD display (1) during the weighing process. The measurement is complete when the display flashes twice.

## Enter or change personal data

**You must enter your personal data when you measure your body weight, BMI, recommended calorie intake and body fat, fluid, bone mass and muscle percentage for the first time.**

### To enter or change personal data:

1. Press **SET** (3) button. Display will show user information.
2. Press **SET**. Memory for user profil (22) flashes in the display (1).
3. Use **UP**-button (2) or **DOWN**-button (4) to choose user profile number (22).
4. Confirm with **SET** (3).

You can then enter in values such as, gender, man (23) or woman (24), age (26), length with **UP** (2) and **DOWN** (4). Confirm each choice with **SET** (3).

## Measuring

**To get comparable measurement results, you should always weigh yourself, if possible, at the same time of day (preferably in the morning). Always weigh yourself, if possible, under the same conditions: an empty stomach, after going to the toilet and without clothes.**

**Stand as steady and still as possible during the measurement for best results. Otherwise, the measurement may be interrupted, and the procedure must be repeated.**

**Note that stone weight unit values are not represented as decimal values. The display is divided into two segments with a colon. Only integers are shown in the st weight unit to the left of the column. The remaining weight is displayed in the weight unit lb to the right of the column.** Example: 6: 8.2 st:lb is shown on the display. This corresponds to 6 stone and 8.2 lbs.

### To perform the measurement, proceed as follows:

1. Press **SET** (3) button.
2. Press **UP** (2) button or **DOWN** button (4) to choose user profile from memory (22).
3. Stand on the scale barefoot (5) and make sure you are standing on the 2 electrodes (6). The person memory indicator flashes (22) a few times before the measurement of body weight begins.
4. After the current weight flashes in the display 2 times, the number 0 will move in the display from left to right and input of other values will begin. Remain on both electrodes (6) until the measurement is complete.

After an approved measurement, the following measurement values will be displayed 2 times in succession:

- Body weight
- Body fat in percent (15) obesity (18), overweight (19), normal weight (20) underweight (21)
- Water in percent (16)
- Bone mass (26)
- Recommended calorie intake per day (17)
- Muscle mass in percent (25)
- BMI (11)
- Weight.

## Switching off

The scale has an automatic switch-off function. After displaying the measured value, the scale will automatically switch off.

## Maintenance

### Change the batteries

**Type:** 1,5V alkaline **Quantity:** 2xAAA

1. Turn the scale upside down. Remove the battery cover.
2. Replace with new AA batteries.  
Make sure the +/- poles are in the correct direction.
3. Put the battery cover back in place.

### Control routines

Own checks can be made with known weights but are not necessary.

Contact your dealer if questions arise.

### Cleaning / Disinfection

**For all maintenance with the scale, the power supply must be disconnected.**

**The scale may only be cleaned with household cleaners or disinfectants, e.g. 70% isopropanol solution.**

**We recommend a disinfectant that is intended for cleaning by damp cloth the surfaces of the scale. Follow the manufacturer's instructions.**

**Do not use polishing agents or strong cleaning agents such as petrol, as they can damage the plastic surfaces.**

**Do not spray the disinfectant on the scale display (with pressure).**

**The disinfectant must not penetrate the scale.**

### Error codes

If the scale does not work as it should or displays error codes, please contact your dealer.

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*Lidén Weighing AB | Sandbergsvägen 3B | 441 39 Alingsås | 0322-442 00*